

# IPV-220 IRRIGATION PUMP

Automatic irrigation system for pneumatic powered instruments used in small bone surgery



## IPS-220 Pump Sensor

Detects small amounts of airflow in the exhaust to start the pump. Fits simply between the air hose and air supply.



## IPV-220 Irrigation Pump


Its design allows the user to quickly and easily select between 4 function settings and adjust flow level. The irrigation tubing clips neatly along the hose to the handpiece to ensure it is kept out of the way.

- Compatible with HST-100 and EPV-220 instruments.
- Fully disposable irrigation kits.
- Simple 'clip-on' set-up.
- 4 function settings.
- Variable flow rate.
- Removeable pump head for cleaning.
- Works with foot or lever control.

IPS-220 pump sensor. Connects to foot pump or air supply via the MA7 probe.

Irrigation tube connected to a De Soutter Medical HST-100 high speed drill.

# SPECIFICATION AND PART NUMBERS

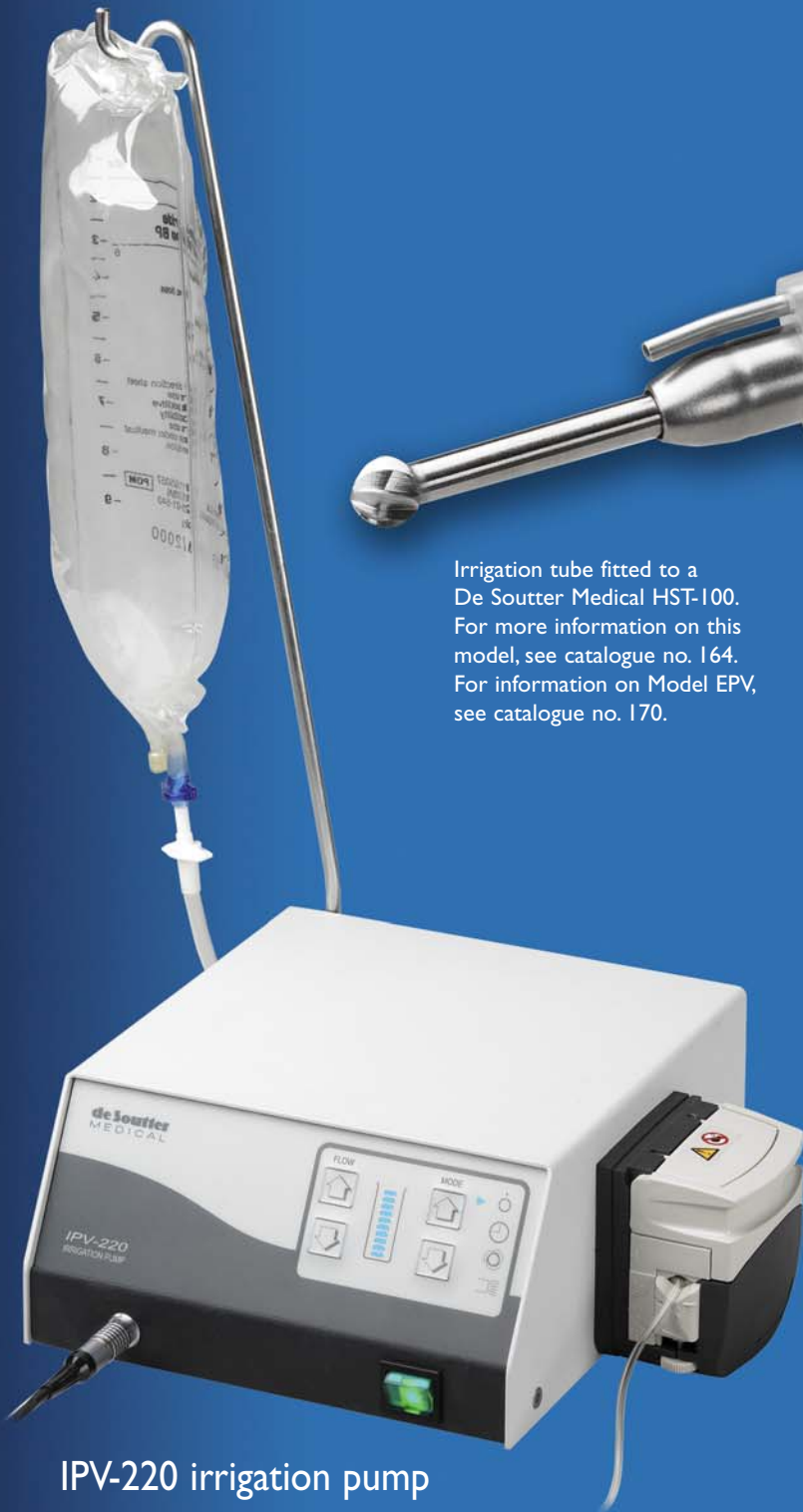

Irrigation tube fitted to a De Soutter Medical HST-100. For more information on this model, see catalogue no. 164. For information on Model EPV, see catalogue no. 170.


## IPV-220 irrigation pump

100 - 240 V                      Order code: 11560  
 Mains power cord sets - 3 m long  
 UK                                      Order code: 8290  
 Euro                                    Order code: 8300  
 US/Canada                        Order code: 8310  
 Australia/N Zealand            Order code: 8320  
 Maximum bag capacity - 1 litre. (bag not supplied)

**De Soutter Medical Limited**  
 Halton Brook Business Park, Weston Road, Aston Clinton, Aylesbury, Bucks HP22 5WF. UK  
 tel: +44 (0) 1296 634 000 fax: +44 (0) 1296 634 033  
 email: info@de-soutter.com web: www.de-soutter.com

**De Soutter Medical Australia Pty Ltd**  
 2/12-14 Apollo Drive, Hallam, Victoria 3803  
 tel: 03 9702 4441 fax: 03 9702 4484  
 email: australia@de-soutter.com web: www.de-soutter.com

## Irrigation pump kits

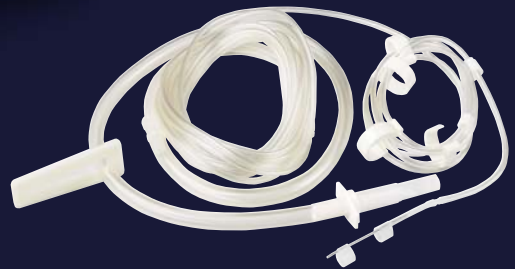

Description	Order code
Irrigation tube set - bur	T10-100
Irrigation tube set - HST-100	T10-102

## Wash Assist hose assemblies



Instrument type	Hose length and order codes		
	3m	4m	5m
EPV, HST	632333	632343	632353

## IPS-220 pump sensor



Order code 11570

**de Soutter**  
 MEDICAL